




novastep®



Forefoot  
Midfoot

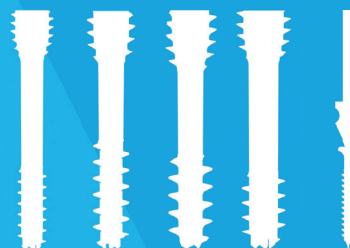

nexis®

## Headless & Snap-off screws


Ø2.3mm | 2.9mm | 3.5mm | 4.0mm | 2.0mm Snap-off

**The Nexis system provides a complete and versatile portfolio of conveniently packaged, easily deployed, low-profile screws designed to address a broad spectrum of forefoot and midfoot indications.**

- Self-Drilling / Self-Tapping
- Reverse cutting flutes
- Torx recess headless screws
- Wire driver compatible snap-off screws
- Self-penetrating conical head
- Simplified color-coded instrumentation
- Sterile packaging
- cleanSTART™ Dispenser

novastep-us.com

# nexis® Headless & Snap-off screws

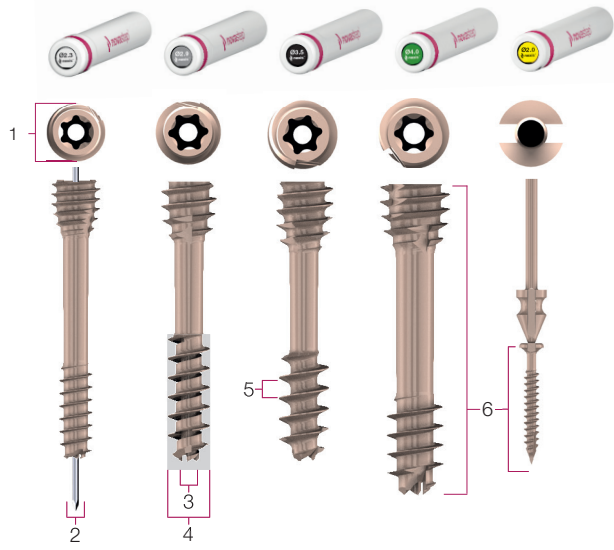


1. Screwdriver Tip	T7 – Torx	T8 – Torx	T8 – Torx	T10 – Torx	Weil†
2. Guide Wire	Ø0.9	Ø0.9	Ø0.9	Ø1.4	N/A
3. Core Diameter	Ø1.75	Ø1.85	Ø1.95	Ø2.65	Ø1.25
4. Thread Diameter	Ø2.3	Ø2.9	Ø3.5	Ø4.0	Ø2.0
5. Thread Pitch	0.9mm	1.35mm	1.45mm	1.35mm	0.85mm
6. Length	10 – 30mm*	10 – 34mm*	14 – 24mm*	18 – 60mm**	11-16mm
7. Drill Bit	Ø1.75††	Ø1.75††	Ø1.75††	Ø2.7††	N/A

\*2mm increments      \*\*2mm increment up to 50 then 5mm increments      †Wire driver compatible      ††Single - use

## Technical Characteristics

Anodized titanium alloy (TA6V). Self-Drilling / Self-Tapping screws with a triple sharpened tip. Reverse cutting flutes aiding the surgeon in simplified screw extraction. Torx recess protects the screw head and provides optimal torque. Self-penetrating conical head design allows atraumatic penetration.

## Quick Tube Deployment

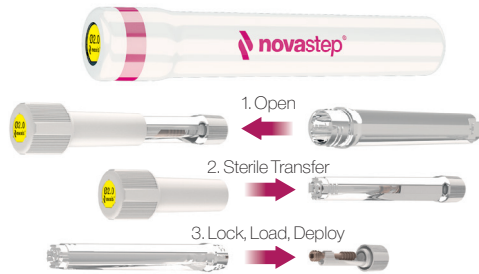


## cleanSTART™ Technology

All implants are packaged sterile and deployed by means of the cleanSTART™ Dispenser. The Nexis screws are ready for immediate use; hospital inventory is easily managed and reliably tracked. Tubes are color coded by screw diameter. Wide choice of lengths ensures a complete portfolio of user friendly implants.



## Simplified Instrumentation

Color coded drivers permit at-a-glance instrument identification. Self-retaining screwdriver tip facilitates manual installation.

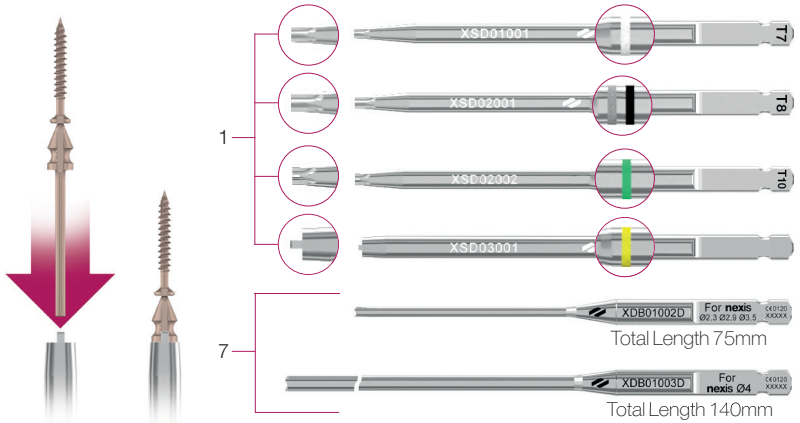

Distributed by:  
Novastep Inc.  
190 Industrial Road, Suite 2  
Wrentham, MA 02093

Manufactured by:  
Novastep SAS  
Espace Performance Alphasis, Bât C1 – C2  
35769 Saint Grégoire, France

For further information call: 877.287.0795

Caution:  
Federal (USA) Law restricts this device for sale to or on the order of a physician. Please refer to package insert for instructions, warnings, contraindications, and potential adverse effects.